

medartis®

PRECISION IN FIXATION

PRODUCT INFORMATION

# Transbuccal Set 2.0/2.3/2.5

MODUS®  
Mandible




# Transbuccal Trocar and Drill Guide

## The MODUS Transbuccal Set

- Compatible with all MODUS sets
- Modular design and time-saving quick release fasteners
- Intra and extraoral cheek retractors

The transbuccal approach to drilling and screw insertion offers a safe osteosynthetic treatment of posterior mandible fractures. The new MODUS transbuccal set features unique flexibly stored drill guides. Even when extracting drilling inserts, the surgeon can apply pressure on the implant and secure its position with the help of the handle. The convenient quick release fasteners and the configuration possibilities with intra or extraoral cheek retractor facilitate the surgeon's task. This new transbuccal set is compatible with all MODUS sets for the treatment of the mandible:

- MODUS 2.0
- MODUS TriLock 2.0/2.3/2.5
- MODUS Trauma 2.5

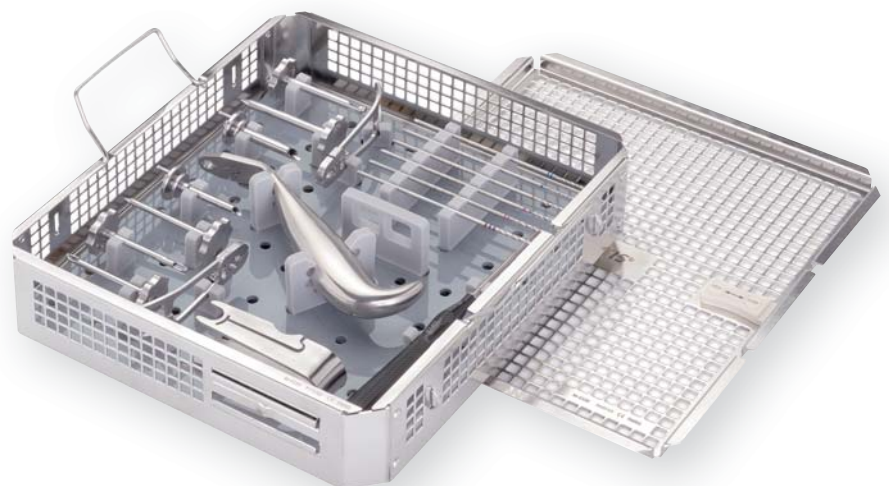

## CONFIGURATION POSSIBILITIES

Protection sleeve Ø 4.5 mm with  
intraoral or extraoral cheek retractor  
for

- MODUS 2.0 screwdriver  
with self-holding blade
- MODUS TriLock 2.0/2.3/2.5  
screwdriver



Protection sleeve Ø 6 mm with  
intraoral or extraoral cheek retractor  
for

- MODUS 2.0 screwdriver  
with tension sleeve
- MODUS 2.5 screwdriver  
with self-holding blade



# Ordering Information

## Handle for Transbuccal Trocar and Drill Guide



Art. No.	Piece/Pkg
M-2606	1

## Cheek Retractors



M-2614



M-2615



M-2616

Art. No.	∅	Description	Length	Piece/Pkg
M-2614		extraoral	166 mm	1
M-2615*	4.5 mm	intraoral	83 mm	1
M-2616*	6.0 mm	intraoral	83 mm	1

\* For fixation to be used with screwdriver M-2598

## Screwdriver



Art. No.	Interface	Description	Length	Piece/Pkg
M-2598	 HD6	self-holding	175 mm	1

## Protection Sleeves



M-2620



M-2621

Art. No.	∅	Description	Length	Piece/Pkg
M-2620	4.5 mm	for M-2606	48 mm	1
M-2621	6.0 mm	for M-2606	48 mm	1

## Trochars



M-2626



M-2627

Art. No.	Description	Length	Piece/Pkg
M-2626	for M-2620	60 mm	1
M-2627	for M-2621	62 mm	1

## Drill Guides



M-2630



M-2631



M-2632

Art. No.	System Size	Description	Length	Piece/Pkg
M-2630	2.0 – 2.5	centric	54 mm	1
M-2631	2.5	eccentric	54 mm	1
M-2632	TriLock 2.0/2.3/2.5		54 mm	1

Twist Drills for Screws Ø 2.0 mm and Cancellous Screws Ø 2.3

Ø 1.5 mm



M-3353



M-3363

Ø 1.6 mm



M-3373



M-3383

Art. No.	Ø Twist Drill	Stop	Length	Drill Shaft End	Piece/Pkg
M-3353	1.5	25 mm	112 mm	Dental	1
M-3363	1.5	25 mm	112 mm	Stryker J-Latch	1
M-3373	1.6	25 mm	112 mm	Dental	1
M-3383	1.6	25 mm	112 mm	Stryker J-Latch	1

Twist Drills for Screws Ø 2.3 mm

Ø 1.9 mm



M-3358



M-3368

Art. No.	Ø Twist Drill	Stop	Length	Drill Shaft End	Piece/Pkg
M-3358	1.9	25 mm	112 mm	Dental	1
M-3368	1.9	25 mm	112 mm	Stryker J-Latch	1

## Twist Drills for Screws Ø 2.5 mm (and for Gliding Hole Ø 2.0 mm)

### Ø 2.0 mm



M-3378



M-3388

Art. No.	Ø Twist Drill	Stop	Length	Drill Shaft End	Piece/Pkg
M-3378	2.0	25 mm	112 mm	Dental	1
M-3388	2.0	25 mm	112 mm	Stryker J-Latch	1

## Twist Drills for Screws Ø 2.8 mm

### Ø 2.3 mm



M-3355



M-3365

Art. No.	Ø Twist Drill	Stop	Length	Drill Shaft End	Piece/Pkg
M-3355	2.3	25 mm	112 mm	Dental	1
M-3365	2.3	25 mm	112 mm	Stryker J-Latch	1

## Twist Drills for Gliding Hole Ø 2.5 mm

### Ø 2.5 mm




M-3375



M-3385


Art. No.	Ø Twist Drill	Stop	Length	Drill Shaft End	Piece/Pkg
M-3375	2.5	25 mm	112 mm	Dental	1
M-3385	2.5	25 mm	112 mm	Stryker J-Latch	1

### Instrument Tray for M-6320

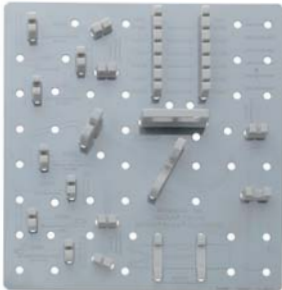

Art. No.	Description	Piece/Pkg
M-6380		1

### Coding Label

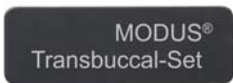


Art. No.	Description	Size	Piece/Pkg
M-6537	MODUS® Transbuccal Set	20 x 60 mm	1

### Instrument Case



M-6320

### Lid for Instrument Case

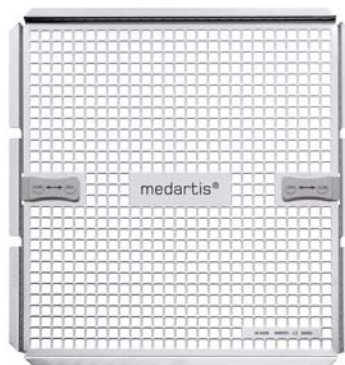

M-6330

Art. No.	Description	Size	Piece/Pkg
M-6320	Stahl	231 x 244 x 54 mm	1
M-6330	Stahl	228 x 243 x 10 mm	1



# System Overview

		SYSTEMS				
		MODUS 2.0 (self-holding)	MODUS 2.0 (with Tension Sleeve)	MODUS TriLock 2.0/2.3/2.5	MODUS Trauma 2.5	MODUS Complete Set
COMPONENTS	Handle	M-2606	M-2606	M-2606	M-2606	M-2606
	Cheek Retractor extraoral	M-2614	M-2614	M-2614	M-2614	M-2614
	Cheek Retractor intraoral	M-2615	M-2616	M-2615	M-2616	M-2615 M-2616
	Tissue Protection Sleeve	M-2620	M-2621	M-2620	M-2621	M-2620 M-2621
	Trocar	M-2626	M-2627	M-2626	M-2627	M-2626 M-2627
	Drill Guide	M-2630	M-2630	M-2632	M-2630 M-2631	M-2630 M-2631 M-2632

TWIN DRILL FOR						
		Dental	M-3353	M-3353	M-3353	
Screws Ø 2.0 < 10 mm	Stryker	M-3363	M-3363	M-3363		M-3363
	Dental	M-3373	M-3373	M-3373		M-3373
Screws Ø 2.0 > 10 mm	Stryker	M-3383	M-3383	M-3383		M-3383
	Dental	M-3378	M-3378	M-3378*		M-3378*
Lag Screws Ø 2.0	Stryker	M-3388	M-3388	M-3388*		M-3388*
	Dental			M-3358		M-3358
Screws Ø 2.3	Stryker			M-3368		M-3368
	Dental			M-3378*	M-3378	M-3378*
Screws Ø 2.5	Stryker			M-3388*	M-3388	M-3388*
	Dental				M-3375	M-3375
Lag Screws Ø 2.5	Stryker				M-3385	M-3385
	Dental				M-3355	M-3355
Screws Ø 2.8	Stryker				M-3365	M-3365

\* Only 1 x in the Set

MANDIBLE-01000001\_v0 / © 01.2011, Medartis AG, Switzerland. All technical data subject to alteration.

#### **HEADQUARTERS**

Medartis AG | Hochbergerstrasse 60E | 4057 Basel/Switzerland  
P +41 61 633 34 34 | F +41 61 633 34 00 | [www.medartis.com](http://www.medartis.com)

#### **SUBSIDIARIES**

Australia | Austria | France | Germany | Mexico | New Zealand | Poland | UK | USA

For detailed information regarding our subsidiaries and distributors, please visit [www.medartis.com](http://www.medartis.com)